



 mitoredlight

MITO **PRO+** | MITO **PRO X** FLOOR STAND

www.mitoredlight.com



ABOUT US

Located in Scottsdale, Arizona, **Mito Red Light**[®] is dedicated to bringing the highest quality red light therapy lights to market at the best value.

Having extensively studied and realized the benefits of red light therapy firsthand, Mito Red Light was brought into existence in order to provide consumers truly affordable, yet highly innovative red light therapy options.

If you have any further questions, please feel free to contact us at:
info@mitoredlight.com

You can also call our phone number
Monday - Friday 9am - 5pm MST
+1 866-861-6486

Mail correspondence can be sent to:
Mito Red Light, Inc
9319 N 94th Way, Suite 400 Scottsdale, AZ 85258

Instructions:

MitoPRO+ / MitoPRO X Floor Stand



Step 1: Unbox your stand and unpackage. You will find:

- 1 stand
- 4 wheel castors
- 2 panel locking screws




Step 2: Thread the provided casters into the base of the floor stand and hand tighten.

Wheels are optional, so skip this step if you prefer not to use them.

Instructions:

MitoPRO+ / MitoPRO X Floor Stand

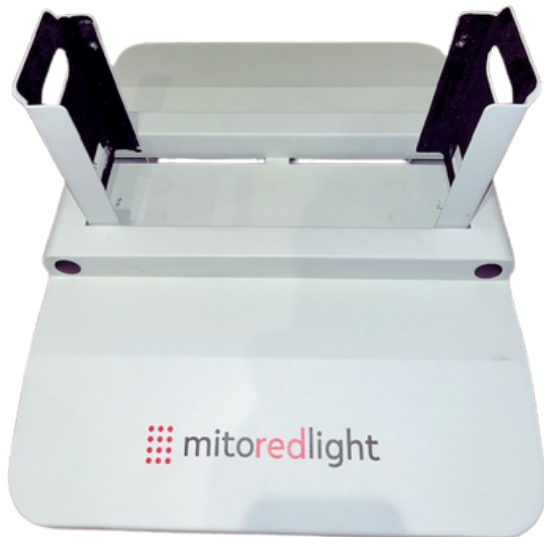

Step 3: Flip the stand onto its casters and flip up the two side braces.



Step 4: Add the two panel locking screws.

Instructions: Panel Placement



Step 1: You will need to remove the lower two rubber feet from the back of your unit using a small Philips head screwdriver before installation into the floor stand.



Step 2: Insert your light into the floor stand and ensure it is seated properly.

Important note: Ensure the panel is fully seated into the bottom of the stand by pushing downwards on the top of the panel.

Instructions: Panel Placement

Step 3: Use the two panel locking screws to secure the panel into place.



Step 4: Add power cord to panel and refer to the panel user manual on how to operate your panel.



Instructions: Wall mount

Your floor stand and panel can be secured to the wall for extra security. Your floor stand package will automatically come with the required equipment.




Step 1: Use the two panel locking screws to secure the panel into place.




Instructions: Wall mount

Step 2: Use the two panel locking screws to secure the panel into place.

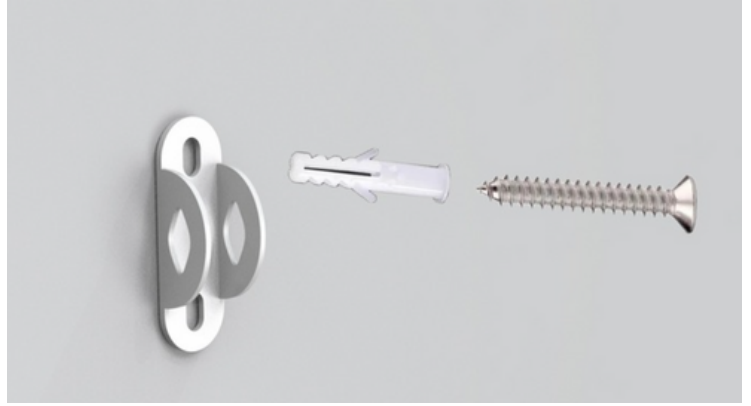

Step 3: Use the two panel locking screws to secure the panel into place.

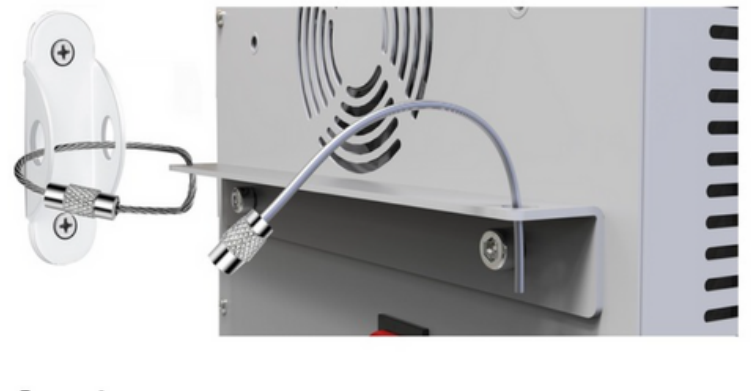

Step 4: Use the two panel locking screws to secure the panel into place.







www.mitoredlight.com

User Manual and Instructional Videos

This user manual and instructional videos can be found here:

<https://mitoredlight.com/pages/user-manual>

Updated January 2025

Follow Us!



[/mitoredlight](https://www.facebook.com/mitoredlight)



[/mitoredlightofficial](https://www.instagram.com/mitoredlightofficial)



[@MitoRedLight](https://www.youtube.com/@MitoRedLight)



[@mitoredlightofficial](https://www.tiktok.com/@mitoredlightofficial)